Aris Hospital

06435 The Aris hospital bed is an efficient and comfortable metal bed. The backrest can be inclined very easily.

Head and foot panel


- The front frame is made of a bent square metal tube. -
- The front frames are firmly bolted to the mattress platform.
- Choice between 2 panel types: a version with a low panel, or a version with a high solid laminate insert panel (*).
- Rounded corners to increase the passive safety of the bed.
- The bottom of the bed panels is equipped with synthetic glides.
- Castors with separate brakes are optionally available.

Mattress platform sections

- One- or two- (*) section metal mattress platform.
- The mattress platform is made of a frame in which rounded steel mattress platform profiles (width: 70 mm) have been welded.
- This construction can be easily cleaned and guarantees a high stiffness and an optimal mattress ventilation.
- The backrest is made of a separate hinged metal frame.
- The backrest can be inclined up to +/- 70° by means of a gas spring.
- The gas spring is activated by a lever that is within reach on both sides of the bed.
- The 2 bed corners at the head end are equipped with cases to install an IV rod or a lifting pole (only available with a two-section mattress platform).



- Width: 90 cm. Length: 207 cm
- Height of the mattress platform: \pm 43 cm
- Safe working load: 200 kg
- Weight (without accessories): 43 kg
- Mattress platform dimensions: 201 x 90 cm (std version)
- Mattress dimensions: 195 x 85 cm Thickness: 12 to 15 cm

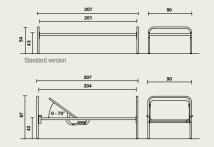

Materials and finishing

- Metal: epoxy-coating RAL 9006 white aluminium.
- Solid laminate.
- Self-lubricating synthetic hinges (= maintenancefree)
- Withstands the common cleaning products.

For more information about the materials and construction: please consult our material lists.

We would to refer to our extensive colour gamut.

Options

- Two-section mattress platform with inclining backrest.
- High bed panels in solid laminate.
- Cases for the Atmosphère sideguards.
- Atmosphère sideguards (only with a solid laminate panel).
- Castors 100 mm with brake (only with a solid laminate panel).

Paneltypes

- Low panel

- High solid laminate insert panel (*)

Sideguards

- Atmosphere (*)

Accessories

We would like to refer to our accessory list.

(*) Option

