# Aris Care

**10010** The Aris Care bed is an efficient and comfortable nursing bed. The optional backrest can be inclined very easily on both sides of the bed.

### Head and foot panel

- The legs are made of solid oven-dried beech.
- The filling panels are made of melaminated chipboard.
- The bed panels are firmly bolted to the mattress platform.

# Mattress platform


- One- or two-section (\*) metal mattress platform.
- The mattress platform is made of a frame in which rounded steel mattress platform profiles of 70 mm have been welded.
- This construction can be easily cleaned and guarantees a high stiffness and an optimal mattress ventilation.
- The backrest is made of a separate hinged metal frame.
- The backrest can be inclined up to +/- 70° by means of a gas spring (only available with a two-section mattress platform).

The legs have a conical shape.Equipped with synthetic glides.

- The gas spring is activated by a lever that is within reach on both sides of the bed.
- The 2 bed corners at the head end are equipped with cases for an IV rod or a lifting pole (only available with a two-section mattress platform).
- The mattress platform sections are equipped with openings for patient restraint straps.

#### Dimensions

- Width: 94 cm. Length: 210 cm
- Height of the mattress platform:  $\pm$  41 cm
- Matress platform dimensions: 204 x 90 cm
  Mattress dimensions: 195 x 90 cm
- thickness: 12 to 15 cm

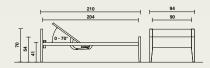

# Materials and finishing

- <u>Wood</u>: Beech (massive, veneer), chipboard (melaminated).
- Metal: epoxy-coating RAL 9006 white aluminium.
- Synthetic material.
- Self-lubricating synthetic hinges (= maintenancefree).
- Withstands the common cleaning products.

For more information about the **materials** and **construction**: please consult our material lists. We would to refer to our extensive **colour gamut**.

#### Options

- Two-section mattress platform with inclining backrest.
- Combination of a low foot panel high head panel
- Wooden surrounds for the mattress platform frame.

#### **Paneltypes**

- Panel Aris Care

#### Accessories

We would like to refer to our accessory list.

(\*) Option

