4ron+

The Aron hospital bed is a height-adjustable multifunctional column bed that provides optimal qualities in the fields of aesthetics, functionality and maintenance.

Carriage

- Frame with 2 anodised aluminium side beams and 2 steel cross beams.
- High/low adjustment by means of 2 telescopic lifting columns.
- Trendelenburg and Reverse-Trendelenburg adjustment up to ±16°.
- Integral S castors Ø 150 mm on precision ball
- bearings, with central brake mechanism.
- Design pedal to activate the central brake mechanism or directional lock.
- The brake pedal is accessible in all bed positions.
- Finished with a synthetic cover cap.
- Open carriage to improve the accessibility under the bed and the ease of cleaning.

Mattress platform

- Mattress platform with 4 sections.
- The backrest is electrically adjustable up to $\pm 70^{\circ}$ with autoregression.
- The backrest is equipped with an inclination indicator.
- The knee rest is electrically adjustable up to $\pm 34^{\circ}$.
- The footrest is manually adjustable for the Fowler
- Electric manipulations of the mattress platform sections are equipped with a clamp safety.
- Manual quick release of the backrest (CPR).
- The perforated, detachable synthetic mattress platform sections can be cleaned very easily.
- Standard equipped with reinforced cutouts for patients restraint straps.
- The backrest and footrest are equipped with integrated mattress holders.

- Telescopic bed extension possible up to 20 cm.
- Equipotential connection.
- 4 bed corners with integrated bumper castors.
- 4 bed corners with cases to install an IV rod or lifting
- The bed panel is made of a metal frame (Vela) or a synthetic frame (* Lauro panel) and a solid laminate
- Aluminium side beams with integrated rail and accessory hooks.
- Built-in cable trough.

Electric bed manipulations

- The electric functions can be activated by means of a hand switch, a mini-nurse box, a nurse box (*) or a satellite control (*).
- The hand switch is equipped with a preprogrammed sleeping function.
- A (mini-)nurse box makes it possible to switch off several functions selectively.

Sideguards

- Vela
- Lauro (*)
- Jaro (* only with Trix)

Panel types

- Senzo (*)
- Atmosphère (*)
- Trix (* only with Jaro panel)
- We would like to refer to our


Accessories

accessory list.




(*) Option

HAELVOET NV | L. BEKAERTSTRAAT 8, 8770 INGELMUNSTER, BELGIUM | T +32 (0) 51 48 66 95 | F +32 (0) 51 48 73 19 | INFO@HAELVOET.BE | WWW.HAELVOET.BE Since Haelvoet continuously aims to make technological progress, Haelvoet reserves the right to modify the product.

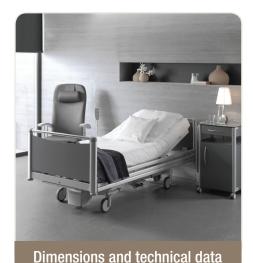

- Width: 100 cm, Length: 211 cm
- Height variation of the mattress platform: ± 41 cm - ± 81 cm (var. 2%)
- Safe working load: 230 kg
- Weight (without options): 140 kg
- Mattress platform dimensions: 200 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 16 cm

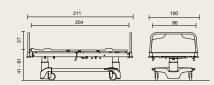

Materials and finishing

- Metal: epoxy-coating RAL 9006 white aluminium, anodised aluminum, chrome (*).
- Synthetic.
- Solid laminate.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the materials and construction: please consult our material lists.

Options

- Removable solid laminate mattress platform sections of 6 mm.
- Footrest extension up to 20 cm.
- Sideguards.
- Fixed head panel.
- Linen holder (not possible with mini-nurse box).
- Nurse box (ACO) with handset + extendible linen holder.
- Nurse box (ACO) with satellite control (FPP) + extendible linen holder.
- Under bed light.
- Double high/low foot pedal with fall protection.
- Electromechanical quick release (CPR).
- Brake buzzer.
- 5th castor.
- Wall spacer.
- Twin castors Ø 150 mm.
- DIN rail (head, foot, head and foot).