

# Artena Hospital

**03873** The Artena hospital bed is a height-adjustable bed that complies with all quality, design and functionality demands of both the staff and the patient.

## Carriage

- Rectangular tube frame.
- Stable high/low movement.
- High-quality Levina castors Ø 125 mm, with separate brakes.
- Double design pedal to activate the central brake mechanism (\*).

## Mattress platform

- Three- or four- (\*) section mattress platform.
- The backrest is electrically adjustable up to  $\pm 70^\circ$ , with auto-regression.
- The knee rest is electrically adjustable up to  $\pm 34^\circ$  (in case of 4 sections).
- The footrest is manually adjustable for the Fowler position (in case of 4 sections).
- The legrest is manually adjustable up to  $\pm 14^\circ$  (in case of 3 sections).
- Electric manipulations of the mattress platform sections are equipped with a clamp safety.
- The mattress platform consists of rounded steel mattress platform profiles.
- Standard equipped with cutouts for patients restraint straps.
- Backrest and footrest are equipped with 2 integrated mattress holders.
- Standard with a telescopic bed extension up to + 20 cm.
- The 4 bed corners are equipped with an integrated bumper.
- The 4 bed corners are equipped with cases to install an IV rod or a lifting pole.
- Choice out of 4 different bed panels (3 with a solid laminate finishing panel and 1 in wood).
- The side beams on both sides are equipped with 3 synthetic accessory hooks.

## Electric bed manipulations

- The electric bed functions can be activated by means of a hand switch.
- The hand switch is equipped with a preprogrammed sleeping function.
- An optional control lock (ACL) (\*) makes it possible to switch off several functions with 1 single button.

### Paneltypes

- Vela (\* lock mechanism)
- Lauro (\*)
- Jaro (\* with Trix sideguards)
- Senzo (\*)

### Sideguards

- Atmosphere (\*)
- Duo-Safe (\*)
- Trix (\* only with Jaro-panel)

### Accessories

We would like to refer to our accessory list.


(\*) Option



## Dimensions and technical data

- Width: 100 cm, Length: 212 cm
- Height variation of the mattress platform<sup>1</sup>:  $\pm 38$  cm -  $\pm 78,5$  cm
- Clearance height<sup>1</sup>: L: 170 mm, R: 135 mm
- Safe working load: 230 kg
- Weight (with options): 125 kg
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 16 cm

<sup>1</sup> With standard model

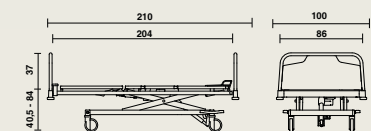

## Materials and finishing

- Metal: epoxy-coating RAL 9006 *white aluminium*, chrome (\*).
- Synthetic.
- Solid laminate.
- Hinge points with self-lubricating synthetic bearings (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials, construction and cleaning products**: please consult our material lists. We would refer to our extensive **colour gamut**.

## Options

- 4-section mattress platform.
- 3- or 4-section mattress platform with removable synthetic sections.
- 3- or 4-section mattress platform with removable solid laminate sections of 6 mm.
- Foot extension up to + 20 cm.
- Panel types: Lauro, Jaro or Senzo.
- Sideguards.
- Total lock device (ACL).
- High/low pedal.
- Equipotential connection.
- Battery.
- Manual quick release of the backrest (CPR).
- Integral S castors Ø 125 mm or Ø 150 mm on precision ball bearings, centrally braked.

