
# Atlas care


The Atlas Care bedside table low model, is very compact and matches perfectly with the Atlas wardrobe gamut.

## Low model (storage cabinet, drawer)

- 2 side panels with a front lath, back, bottom, and a tabletop.
- Made of melaminated chipboard (thickness 18 mm).
- The panels are connected to each other by means of wooden dowels and eccentric cams.
- Storage cabinet with a removable shelf.
- High-quality self-closing hinges, opening up to 135°.
- Metal design handle.
- The removable drawer is made of metal sides, between which a wooden bottom and back have been attached.
- The drawer is installed in the cupboard frame by means of high-quality steel rails with an end stop.
- The metal parts are protected by an epoxy-coating.
- The tabletop has rounded corners and edges and is made of melaminated chipboard of 18 mm.
- 4 adjustable glides.

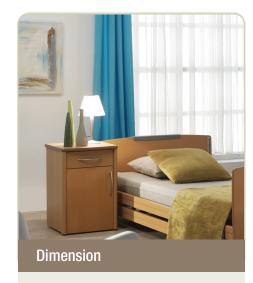

**06241** | 44 x 44 x height 63 cm

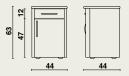

### Materials and finishing

- Wood: chipboard (melaminated).
- Metal: epoxy-coating (clear gray).
- High-quality self-closing hinges.
- Withstands the common cleaning products.

For more information about the materials, construction and maintenance: please consult our material lists.

### **Options**

- Finishing in veneer.
- Cases for accessories on the tabletop.
- Telephone holder.
- 4 hidden castors.

# Colours

We would like to refer to our extensive colour gamut. Sample card + samples available on request.

