

Atlas Hospital

07169 The Atlas Hospital is a compact bedside table that offers numerous storage possibilities thank to its optimal design.

One-sided version (storage cabinet, open compartment, drawer)

- 2 side panels with a front lath, back, bottom, and a tabletop.
- The frame, door, and drawer fronts are made of light grey melaminated chipboard (thickness 18 mm).
- The panels are connected to each other by means of wooden dowels and eccentric cams.
- The front posts have a melamine layer with an aluminium look.
- Storage cabinet with a removable shelf.
- High-quality self-closing hinges, opening up to 135°.
- Metal design handle.
- The removable drawer is made of metal sides, between which a wooden bottom and back have been attached. The drawer is finished with a wooden drawer front.
- The drawer is installed in the cupboard frame by means of high-quality steel rails with an end stop.
- The metal parts are protected by a light grey epoxy-coating.
- The tabletop has rounded corners and edges and is made of a melaminated chipboard of 18 mm.
- 4 hidden castors Ø 50 mm.


Optional dinner tray Livia

- Folding.
- Easy to incline in both directions.
- The dinner tray has an upright edge.
- The folding part is fixed on a strong aluminium guiding profile that has been finished with an anodisation layer.
- The dinner tray is height-adjustable (between 80 and 110 cm) by means of a gas spring.
- Dimensions of the dinner tray: 63 x 40,5 cm.
- The gas spring is activated by means of a foot pedal.
- The dinner tray is made of fire-retardant ABS.
- The bedside table stands on 5 design castors Ø 65 mm.



Dimension

07169 | 49,5 x 49,5 x height 83 cm
with dinner tray Livia | 65 x 49,5 x height 88 cm

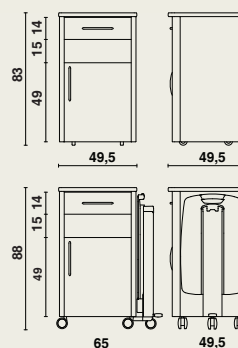

Materials and finishing

- Wood: chipboard (melaminated, light grey).
- Metal: epoxy-coating (clear gray).
- High-quality self-closing hinges.
- Withstands the common cleaning products.

For more information about the materials, construction and maintenance: please consult our material lists.

Options

- Livia dinner tray.
- Gallery.
- Cylinder lock.
- Cases for a telephone holder.
- Telephone holder.
- 2nd drawer in the open compartment.
- Extensible tray.
- 4 design castors Ø 65 mm.
- 4 design castors Ø 65 mm with brake.

