

# Commode chair wood

**00958** This stable commode chair has a double discrete seat and is therefore very suitable for patient rooms.

## Support

- The frame is made of 2 solid beech sides.
- The sides are connected with a solid beech cross beam at the back.
- The armrests are milled ergonomically.
- A high solid beech beam in front hides the bucket.
- The wooden parts are finished with a coat of varnish.
- The legs are equipped with synthetic glides.
- The chair can be optionally equipped with a metal support and 4 castors with brakes  $\varnothing$  75 mm.


## Seat and backrest

- The double seat consists of a tiltable upper part and a tip-up toilet seat underneath.
- The upholstered seat can be tilted by means of 2 rust-proof hinges or shoved backwards by means of 2 stainless steel brackets.
- The varnished plywood toilet seat with oval opening is screwed onto the frame by means of stainless steel hinges.
- A removable synthetic bucket with a lid is installed underneath the tip-up toilet seat. This bucket rests on stainless steel suspension brackets.
- The seat is upholstered with fire-retardant foam (thickness: 20 mm; density: 50 kg/m<sup>3</sup>).
- The upholstered backrest is made of 12 mm ergonomic plywood that is firmly screwed by means of fixing plates.
- The backrest is equipped with fire-retardant foam (thickness: 10 mm; density: 35 kg/m<sup>3</sup>).



## Dimension

**00958** | 50 x 46 x h 82,5 cm

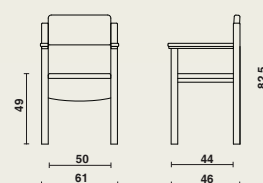

## Materials and finishing

- Wood: beech (Laminated, massive), plywood, MDF.
- Metal: epoxy-coating (light gray).
- Fire-retardant foam.
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials**, **constructions** and **maintenance**: consult our materialfiles.

## Options

- Tip-up solid laminate toilet seat.
- Support with 4 castors  $\varnothing$  75 mm.

## Colours

We would like to refer to our extensive colour gamut.  
Sample card + samples available on request.

