

Domus Classic Flex

08901 The central idea behind the Domus relax chair is to offer optimal seating comfort in peace and quiet.

The Classic mechanism can be easily adjusted by pushing yourself off. Lever manipulations belong to the past.

Support

- The rounded sides consist of 2 OSB boards of 18 mm, and are firmly connected with a metal structure.
- The sides are upholstered with fire-retardant foam.
- The upholstery is finished with imitation leather and decorative stitching.
- Comfortable rounded armrests of ± 5 cm.
- All edges and corners have been rounded.
- No visible fastening elements.
- Design castors $\varnothing 65$ mm with brake.
- Optionally possible on solid beech legs, with aid castors at the back.

Seat

- Removable seat in melaminated MDF, thickness: 12 mm.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 45 cm.
- Upholstered with fire-retardant foam. Two different density zones guarantee optimal seating comfort.
- The upholstery is finished with imitation leather and decorative stitching.

Backrest

- The backrest consists of an ergonomic pre-shaped plywood shell, thickness 12 mm.
- Upholstered with fire-retardant foam. This foam has 3 different density zones to guarantee optimal seating comfort.
- The upholstery is finished with imitation leather and decorative stitching.
- Equipped with an adjustable and removable head/neck cushion.

Legrest

- Plywood, thickness: 10 mm.
- Upholstered with fire-retardant foam and finished with imitation leather.

Adjustable positions

- 3 positions: seating position, relax position and relax position.
- The relax mechanism first raises the legs, after which it is possible to further incline the entire chair.
- Easy to adjust by pushing yourself off.
- Adjustable without lever manipulation.
- Optimal seating comfort in each position, thanks to the constant angle between the backrest and seat during the inclination.
- Optimal leg support, thanks to the telescopic legrest.
- In the seating position, the legrest is turned inward to make it easier to get up.

Colours


We would like to refer to our extensive colour gamut.
Sample card + samples available on request.



Dimension

08901 | Domus Classic Flex:

- L 62 x D 73 x H 115 cm
- Seat height: 45 cm
- Seat depth: 51 cm
- Seat width: 52 cm
- Height armrests: 63 cm

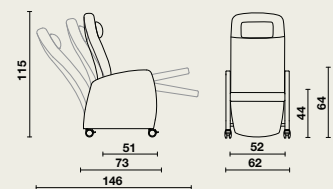

Materials and finishing

- Wood: Solid beech, Plywood, OSB, MDF.
- Metal: epoxy coating (light grey).
- Fire-retardant foam.
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials**, **constructions** and **maintenance**: consult our materialfiles.

Options

- Wooden armrests.
- Dinner tray attachment.
- Abduction block combined with the dinner tray.
- Cushion with visco-elastic foam and a bi-elastic cover.
- Anatomical reduction cushion.
- Lateral headrest.
- Armrests with removable upholstery.
- Case for IV rod.
- IV rod.
- Wooden legs and aid castors at the back.
- Design castors $\varnothing 100$ mm with brake.
- Design castors $\varnothing 100$ mm with brake and 2 directional castors.
- Footrest.

