

Olympia Care

08888 The Olympia nursing bed: well-ried, durable and safe. Thanks to several bed panel types, an extensive colour gamut and the possibility to hide the castors from view, the Olympia Care allows you to give a personal touch to the living room.

Carriage

- Very stable high/low adjustment by means of supporting arms.
- Trendelenburg (*) and Reverse Trendelenburg adjustment up to $\pm 17^\circ$.
- High-quality design castors \varnothing 125 mm, with separate brakes.
- High-quality design castors Integral S \varnothing 125 or 150 mm with central brake (*).
- Design pedal to activate the central brake mechanism (*).
- Hidden central braked castors (\varnothing 50 mm) are available by means of synthetic cover caps (*).
- The carriage has no actuators, which facilitates its accessibility and its possibility to clean.

Mattress platform

- Three- or four- (*) section mattress platform.
- The backrest is electrically adjustable up to $\pm 70^\circ$, with auto-regression.
- The knee rest is electrically adjustable up to $\pm 34^\circ$ (in case of 4 sections).
- The footrest is manually adjustable to get the Fowler position (in case of 4 sections).
- The legrest is manually adjustable up to $\pm 14^\circ$ (in case of 3 sections).
- The electric manipulation of the mattress platform sections is equipped with a clamp safety.
- The mattress platform consists of rounded steel mattress platform profiles.
- The mattress platform sections are equipped with openings for patient restraint straps.
- Telescopic bed extension possible up to 20 cm (not with Pollux sideguards).
- The head end is equipped with 2 cases to install an IV rod or a lifting pole.
- No visible fastening bolts in the foot panel.
- The bed panels are made of solid oven-dried beech (*).
- Rounded corners to increase the passive safety.

Electric bed manipulations

- The electric functions can be activated by means of a hand switch or a *mini-nurse box* (*).
- The hand switch is equipped with a preprogrammed sleeping function.
- A *mini-nurse box* (*) makes it possible to switch off several functions selectively, whereas a lock box (*) switches these functions off by means of 1 single button.
- The optional hand switch (*) can switch off several functions selectively by means of a magnetic key.

Paneltypes

- Midea (standard)
- Valentino (*)
- Puro (*)
- Felina (*)

Sideguards

- Pollux (*)
- Trix (*)
- Atmosphere (*)

Accessories

We would like to refer to our accessory list.



(*) Option




HAELVOET NV | L. BEKAERTSTRAAT 8, 8770 INGELMUNSTER, BELGIUM | T +32 (0) 51 48 66 95 | F +32 (0) 51 48 73 19 | INFO@HAELVOET.BE | WWW.HAELVOET.BE
Since Haelvoet continuously aims to make technological progress, Haelvoet reserves the right to modify the product.



Dimensions and technical data

- Width: 100,5 cm, Length: 208 cm
- Height variation of the mattress platform:
 $\pm 35,5^*$ cm - $\pm 78,5$ cm (with castors \varnothing 125 mm)
 $\pm 34,5^*$ cm - $\pm 77,5$ cm (with hidden castors)
- Clearance height of the carriage:
L: 170 mm, R: 135 mm (with castors \varnothing 125 mm)
L: 160 mm, R: 125 mm (with hidden castors)
- Safe working load: 230 kg
- Weight (without options): 140 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm
(* with Trendelenburg + 2,5 cm)

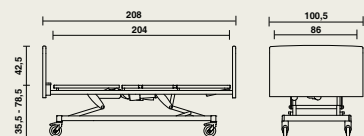

Materials and finishing

- Wood: Beech (massive, veneer), chipboard (melaminated).
- Metal: epoxy-coating RAL 9006 white aluminium.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials and construction**: please consult our material lists. We would refer to our extensive colour gamut.

Options

- 4-section mattress platform.
- Removable synthetic mattress platform sections.
- Foot extension possible up to + 20 cm.
- Panel types: Puro, Felina or Valentino.
- Bumpers and/or wall spacer.
- Sideguards.
- Total lock device (ACL).
- Mini-nurse box (ACC).
- Hand switch with selective locking system (incl. Trend).
- Illuminated hand switch and/or nightlight under the bed.
- Wooden surrounds for carriage and/or mattress platform frame
- Accessory hooks.
- Battery.
- Manual quick release of the backrest (CPR).
- Linen holder.