# Jlympia XLow

10785 The Olympia Care XLow bed is a perfect symbiosis of safety and comfort. A low bed position of  $\pm$  25 cm reduces the risk of falling to an absolute minimum, whereas the ergonomic working height of  $\pm$  80 cm guarantees efficient care.

# Carriage

- Highly stable high/low adjustment by means of supporting arms.
- Trendelenburg (\*) and reverse Trendelenburg adjustment up to  $\pm$  16°.
- Hidden castors Ø 50 mm with synthetic cover caps, braked 2 per 2.
- Design pedal to activate the brake mechanism.
- The carriage has no actuators, which facilitates its accessibility and possibility to clean.
- The use of a patient lift is perfectly possible.
- The bed can be moved in every position, except in the safety position (= lowest position).

### **Mattress platform**

- Three- or four- (\*) section mattress platform.
- The backrest is electrically adjustable up to  $\pm 70^{\circ}$ , with auto-regression.
- The knee rest is electrically adjustable up to  $\pm 34^{\circ}$  (in case of 4 sections)
- The footrest is manually adjustable to obtain the Fowler position (in case of 4 sections).
- The legrest is manually adjustable up to  $\pm 14^{\circ}$  (in case of 3 sections).
- The electric manipulation of the mattress platform sections is equipped with a clamp safety.
- The mattress platform consists of rounded steel mattress platform profiles.

- The mattress platform is equipped with openings for patient restraint straps.
- The backrest and footrest are equipped with 2 synthetic mattress holders.
- Telescopic bed extension possible up to 20 cm (not with Pollux sideguards).
- The head end is equipped with 2 cases to install an IV rod or a lifting pole.
- No visible fastening bolts in the foot panel.
- The bed panels are made of solid oven-dried beech (\*).
- Round corners to increase the passive safety.

#### Electric bed manipulations

- The electric functions can be activated by means of a hand switch.
- The hand switch is equipped with a preprogrammed sleeping function.

Sideguards

Pollux (\*)

Trix (\*)

- The hand switch (\*) can switch off several functions selectively by means of a magnetic key.

We would like to refer to our

### **Paneltypes**

- Midea (standard)
- Valentino (\*)
- Puro (\*)
- Felina (\*)

#### Accessories

accessory list.


(\*) Option





# Dimensions and technical data

- Width: 100,5 cm, Length: 208 cm
- Height variation of the mattress platform: ± 25 cm - + 80 cm
- Clearance height of the carriage: 150 mm (with a mattress platform height of 42 cm)
- Safe working load: 230 kg
- Weight (without options): 140 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm

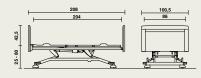

## Materials and finishing

- Wood: Beech (massive, veneer), chipboard (melaminated).
- Metal: epoxy-coating RAL 9006 alu.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products and disinfectants.

For more information about the materials and construction: please ask for our material lists. We would like to refer to our extensive colour gamut.

#### **Options**

- 4-section steel mattress platform.
- Mattress platform with removable synthetic sections.
- Footrest extension up to 20 cm.
- Panel types Puro, Felina or Valentino.
- Bumpers (not possible with Midea)
- Wall spacer.
- Sidequards.
- Illuminated hand switch.
- Hand switch with selective locking system (incl. Trendelenburg).
- Reading lamp.
- Battery.
- Central brake castors Ø 125 mm activated by a brake bracket.
- Central brake castors Ø 100 mm activated by a brake bracket.