The Olympia hospital bed combines simplicity, elegance and functionality in an utmost ergonomic hospital bed.

08950 | Olympia Hospital II

Electric protection class II

The bed is standard equipped with a double-insulated control box, without earthing, equipotential connection, quick release and battery.

Carriage

- Very stable high/low adjustment by means of a support arm system.
- Trendelenburg (*) and Reverse Trendelenburg adjustment up to $\pm 17^{\circ}$.
- Individually braked design castors Ø 125 mm.
- Integral S castors Ø 125 mm or Ø 150 mm on

precision ball bearings, centrally braked (*).

- Synthetic design pedal to activate the central brake mechanism (*).
- The open carriage improves the accessibility under the bed and the ease of cleaning.

Mattress platform

- Matress platform with 3 or 4 (*) sections.
- The backrest is electrically adjustable up to $\pm 70^{\circ}$ with autoregression.
- The knee rest is electrically adjustable up to $\pm 34^{\circ}$ (in case of 4 sections).
- Legrest or footrest can be manually adjusted.
- The electric manipulation of the mattress platform section is equipped with a clamp safety.
- Mattress platform with rounded steel profiles.
- Standard equipped with cutouts for patients restraint
- Backrest and footrest are equipped with 2 integrated

- mattress holders.
- Standard with a telescopic bed extension up to +
- The 4 bed corners are equipped with integrated bumpers.
- The 4 bed corners are equipped with cases to install an IV rod or a lifting pole.
- Choice out of 3 different bed panels with a solid laminate filling panel.
- The side beams on both sides are equipped with 3 synthetic accessory hooks.

Electric bed manipulations

- The electric functions can be activated by means of a hand switch, a mini nurse box (*) or a nurse box (*).
- The hand switch is equipped with a preprogrammed sleeping function.
- A mini nurse box (*) or nurse box (*) makes it possible to switch off several functions selectively.
- The control lock (*) blocks all electric functions of the hand switch by means of one single button.

Panel types

- Vela
- Lauro (*)
- Jaro (*)
- Senzo (*)

Sideguards

- Atmosphere (*)
- Duo-Safe (*)
- Trix (only with Jaro-panel) (*)

Accessoiries


We refer to our extensive accessory list





Dimensions and technical data

- Width: 100 cm, Length: 212 cm
- Height variation of the mattress platform1: ± 38 cm - ± 78,5 cm
- Clearance height1: L: 170 mm, R: 135 mm
- Safe working load: 230 kg
- Weight (with options): 125 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm
- ¹ With standard model

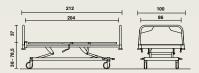

Materials and finishing

- Metal.
- Solid laminate.
- Epoxy-coating RAL 9006 white aluminium or
- Hinge points with self-lubricating synthetic bearings (= maintenance-free).
- Withstands the common cleaning products.

For more information about the materials and construction: please consult our material lists. We would to refer to our extensive colour gamut.

Options

- 4-section matress platform.
- Removable synthetic mattress platform sections.
- Removable solid laminate mattress platform sections of 6 mm.
- Foot extension possible up to + 20 cm.
- Panel types: Lauro, Jaro or Senzo.
- Sidequards.
- Equipotential connection.
- Control lock (ACL) and Mini-nurse box (ACO).
- Nurse box (ACO) + extendible linen holder.
- Illuminated hand switch.
- Nightlight under the bed.
- Battery with alarm signal.
- Trendelenburg adjustment
- Manual quick release of the backrest.
- 5th castor.
- High/low pedal.
- Choice of different type of castors + central brake.



(*) Option