# Quadra Care


The Quadra Care bedside table is very compact and has a contemporary and rigid design. The extra low bedside table matches perfectly with the XLow bed range.

### Low model with 2 drawers

- 2 side panels with a back, bottom, and a tabletop.
- Made of melaminated chipboard (thickness 19 mm).
- The side panels, the top of the cabinet and the bottom The drawer is installed in the cupboard frame by are connected to each other by means of wooden dowels. Then cold glue is used to create a strong cupboard frame.
- The tabletop is connected to the cabinet by means of eccentric cams.
- Metal design handle.
- The removable drawer is made of metal sides,


between which a wooden bottom and back have been attached.


- means of high-quality steel rails with an end stop.
- The tabletop is made of a double sided melaminated chipboard of 19 mm.
- 4 hidden castors Ø 50 mm.

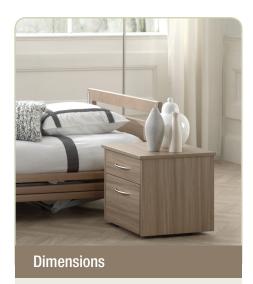

Quadra XLow hidden castors | 48 x 47 x height 46

Quadra XLow castors Ø 65 mm | 48 x 47 x height

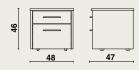

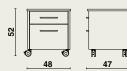

## Materials and finishing

- Wood: chipboard (melaminated).
- Withstands the common cleaning products.

For more information about the materials, construction and maintenance: please consult our material lists.

# **Options**

- Cylinder lock.
- 4 castors Ø 65 mm.
- 4 castors Ø 65 mm with brake.

### Colours

We would like to refer to our extensive colour gamut. Sample card + samples available on request.

