# Siena Hospital

The Siena hospital bed offers you a very functional and efficient mattress platform. Thanks to its flexible finishing choice, this bed can be perfectly combined with other hospital beds from our gamut.

# Carriage

- Legs at a fixed height or adjustable in 3 positions (\*). The carriage can be equipped with a central brake
- The bed can be equipped with 4 design castors Ø 100 mm, with separate brakes (\*).
- mechanism and 4 design castors Ø 125 mm (not height-adjustable) (\*).

# **Mattress platform**

- Three- or four- (\*) section mattress platform, manually adjustable.
- Available with a mattress platform electrically adjustable (\*).
- The backrest is adjustable up to ±70°, with autoregression by means of a gas spring or an electric
- The knee rest is adjustable up to  $\pm 34^{\circ}$  by means of a gas spring (in case of 4 sections).
- The legrest or footrest can be manually adjusted.
- The mattress platform consists of rounded steel mattress platform profiles.
- Standard equipped with cutouts for patients restraint straps

- The backrest and footrest are equipped with 2 synthetic mattress holders.
- Standard equipped with a telescopic bed extension up to + 20 cm
- The 4 bed corners are equipped with cases to install an IV rod or a lifting pole.
- Choice out of 4 different panel types (3 with a solid laminate filling panel and 1 in wood).
- The side beams on both sides are equipped with 3 synthetic accessory hooks.

# **Bed manipulations**

- The backrest and knee rest are supported by a gas spring and can be adjusted by means of a lever on both bed sides.
- The electric back- and knee rest is activated by a hand switch (\*).
- The electric finishing can be equipped with a lock box, making it possible to switch off the electric functions with 1 single button (\*).

### **Sidequards Paneltypes**

Atmosphere (\*)

Trix (\* only with Jaro)

- Vela (standard)
- Lauro (\*)
- Jaro (\* and only with Trix)
- Senzo (\*)

(\*) Option


## Accessories

We would like to refer to our accessory list.



# Dimensions and technical data

- Width: 100 cm, Length: 208 cm
- Height of the mattress platform: ± 46 cm, adjustable up to  $\pm$  55 cm (in 3 positions) (\*).
- Safe working load: 180 kg
- Weight (without options): 65 kg
- Matress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm

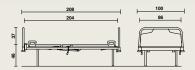

# Materials and finishing

- Metal: epoxy-coating RAL 9006 white aluminium, chrome (\*).
- Solid laminate.
- Self-lubricating synthetic hinges (= maintenance-
- Withstands the common cleaning products.

For more information about the materials and construction: please consult our material lists. We would to refer to our extensive colour gamut.

# **Options**

- 4-section steel matress platform.
- 3- or 4- section mattress platform with removable solid laminate sections of 6 mm.
- Foot extension up to + 20 cm.
- Panel types: Lauro, Jaro and Senzo.
- Wall spacer.
- Sidequards.
- Total lock device (ACL).
- Electrical hand switch.
- Trendelenburg function. High/low cariage Vico.
- Removable vacuum-formed synthetic cover cap.
- Legs adjustable in 3 positions.
- 4 design castors Ø 100 mm, with separate brakes.
- 4 design castors Ø 125 mm with central brake (not height-adjustable).

