# Siena care

The Siena nursing bed offers you a very functional and efficient mattress platform. Thanks to its flexible finishing choice, this bed can be perfectly combined with other nursing beds from our gamut.

# Carriage

- Legs at a fixed height or optionally adjustable in 3
- The bed can be optionally equipped with 4 design castors Ø 100 mm, with separate brakes.
- The carriage can be optionally equipped with a central brake mechanism and 4 design castors Ø 125 mm (not height-adjustable).

# **Mattress platform**

- Three- or four- (\*) section mattress platform, manually adjustable.
- The backrest is adjustable up to  $\pm 70^{\circ}$ , with autoregression by means of a gas spring of electric
- The knee rest is adjustable up to  $\pm$  34° by means of a gas spring (in case of 4 sections).
- Available with a mattress platform electrically
- The legrest or footrest can be manually adjusted.
- The mattress platform consists of rounded steel mattress platform profiles.

- Available with removable solid laminate mattress platform sections of 6 mm (\*).
- The backrest and footrest are equipped with 2 synthetic mattress holders.
- Standard equipped with a telescopic bed extension up to +20 cm.
- The head end is equipped with 2 cases to install an IV rod or lifting pole.
- No visible fastening bolts in the foot panel.
- The bed panels are made of solid oven-dried beech (\*).
- Rounded corners to increase the passive safety.

# **Bed manipulations**

- The backrest and knee rest are supported by a gas spring and can be adjusted by means of a lever on both bed sides.
- The electric back- and knee rest (\*) is activated by a hand switch.
- The electric finishing can be equipped with a lock box (\*), making it possible to switch off the electric functions with 1 single button.

## **Sidequards Paneltypes**

- Midea (standard)
- Valentino (\*)
- Puro (\*)
- Felina (\*)

- Pollux (\*)
- Trix (\*)
- Atmosphere (\*)


# Accessories

We would like to refer to our accessory list.





- Width: 100,5 cm, Length: 208 cm - Height of the mattress platform: ± 46 cm, optionally adjustable up to  $\pm$  55 cm (in 3 positions).
- Safe working weight: 180 kg
- Weight (without options): 65 kg
- Matress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm

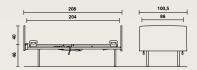

# Materials and finishing

- Wood: Massive beech, melaminated chipboard.
- Metal: epoxy-coating RAL 9006 white aluminium.
- Synthetic material
- Self-lubricating synthetic hinges (= maintenance-
- Withstands the common cleaning products.

For more information about the materials and construction: please consult our material lists. We would to refer to our extensive colour gamut.

# **Options**

- 4-section steel matress platform.
- Mattress platform with removable solid laminate sections.
- Foot extension up to + 20 cm.
- Panel types: Puro, Felina or Valentino.
- Bumpers and/or wall spacer.
- Sideguards.
- Total lock device (ACL).
- Electrical hand switch.
- High/low cariage Vico
- Removable vacuum-formed synthetic cover cap.
- Wooden surrounds of the mattress platform frame.
- Wooden surrounds for the legs.
- Accessory hooks.
- Blanket rack.
- Legs adjustable in 3 positions.
- 4 design castors Ø 100 mm, with separate brakes.
- 4 design castors Ø 125 mm with central brake (not heightadiustable).

