

Talo Classic with castors

12901 The Talo Classic with castors is a comfortable and well-designed relax chair. The Classic mechanism can be easily adjusted by pushing yourself off a little from the chair, and this without using any levers.

Support

- The support is made of 2 solid beech sides, connected to each other by means of a solid beech cross beam.
- The side parts are connected to each other by means of a glued mortise and tenon joint.
- The sides are firmly connected to each other by means of a solid beech cross beam and a metal frame.
- The ergonomic design of the armrests makes it easy to stand up or sit down.
- Comfortable rounded armrests of ± 7 cm.
- All edges and corners have been rounded.
- Four design twin castors \varnothing 100 mm with separate brakes.

Seat

- Removable seat, made of a melaminated MDF board of 12 mm.
- The rounded seat in front stimulates a fluent blood circulation.
- Comfortable seat height of 47 cm.
- Upholstered with fire-retardant foam in 2 different density zones to guarantee an excellent seat comfort.
- The upholstery is finished with imitation leather and decorative stitching.

Backrest

- Made of an ergonomically preformed plywood shell of 12 mm.
- Upholstered with fire-retardant foam in 3 different density zones to guarantee an excellent seat comfort.
- The upholstery is finished with imitation leather and decorative stitching.
- Equipped with an adjustable and removable head/neck cushion.
- Large, ergonomic push bar made of synthetic material at the back.

Legrest

- Made of a plywood board of 10 mm.
- Upholstered with fire-retardant foam and finished with imitation leather.

Classic mechanism

- 3 positions: seat position, reading position, and a relax position.
- The relax mechanism first raises the legs, after which it is possible to further incline the entire chair.
- Easy to adjust by pushing yourself off a little.
- Adjustable without needing a lever.
- Guaranteed optimal seat comfort in every position: when the chair is reclined, the angle between the backrest and seat remains constant.
- Optimal leg support, thanks to the telescopic legrest.
- When the chair is put in the seat position, the legrest is turned inwards so that you can stand up without any problem.


Colours

We would like to refer to our extensive colour gamut.
Sample card + samples available on request.

Dimensions

12901 | W 68 x D 95 x H 124 cm

- Seat height: 47,5 cm
- Seat depth: 51 cm
- Seat width: 53,5 cm
- Height armrests: 65,5 cm

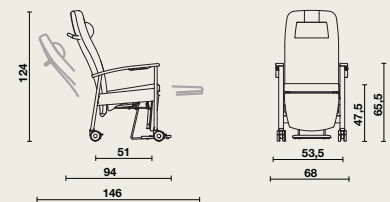

Materials and finishing

- Wood: beech (solid), plywood, MDF.
- Metal: epoxy coating (light grey).
- Fire-retardant foam.
- Upholstery: see collection.
- Resistant to the common cleaning products.

For further information about used **materials, constructions and maintenance**: consult our materialfiles.

Options

- Dinner tray attachment on armrests.
- Seat cushion with visco-elastic foam and a bi-elastic cover.
- Anatomical reduction cushion.
- Lateral headrest.
- Backrest with removable headrest.
- Closed side panels.
- Accessory brackets.
- Restraint bracket.
- Case for IV rod.
- Dinner tray support on backrest.
- Footrest.

