

Vico Care

08750 The Vico nursing bed combines well-trying technology with a homely look. The hydraulic high/low adjustment guarantees maximum stability, flexibility and safety.

Carriage

- Hydraulic high/low adjustment by means of support arms.
- The steel frame is finished with a removable vacuum-formed synthetic cover cap (*).
- Trendelenburg adjustment up to $\pm 12^\circ$ (*).
- Reverse Trendelenburg adjustment up to $\pm 13^\circ$ (*).
- High-quality design castors \varnothing 125 mm, with separate brakes.
- Design pedal to activate the central brake mechanism (*).

Mattress platform

- Three- or four- (*) section mattress platform.
- Adjustable backrest up to $\pm 70^\circ$ with auto-regression by means of a gas spring.
- Adjustable knee rest up to $\pm 34^\circ$ by means of a gas spring (only with 4 sections).
- Legrest or footrest are manually adjustable.
- Equipped with rounded steel mattress platform profiles.
- Available with removable solid laminate mattress platform sections of 6 mm (*).
- Backrest and footrest are equipped with 2 synthetic mattress holders.
- Standard equipped with a telescopic bed extension up to + 20 cm.
- The head of the bed is equipped with 2 cases to install an IV rod or lifting pole.
- No visible fastening bolts on the foot panel.
- Bed panel available in solid oven-dried beech (*).
- Rounded corners to increase the passive safety.

Electric bed manipulations

- High/low adjustment by means of pedals within easy reach on both sides of the bed.
- The Trendelenburg position is activated by means of the pedal at the foot of the bed (*).
- The gas spring-supported backrest and knee rest can be adjusted by means of a lever on both sides of the bed.

Paneltypes

- Midea (standard and only with Pollux)
- Valentino (*)
- Hypnos (*)
- Felina (*)

Sideguards

- Pollux (*)
- Trix (*)
- Atmosphere (*)

Accessories


We would like to refer to our accessory list.



Dimensions and technical data

- Width: 100,5 cm, Length: 208 cm
- Height variation of the mattress platform¹: ± 41 cm - ± 80 cm
- Safe working weight: 180 kg
- Weight (without accessories): 115 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm

¹ For a standard model

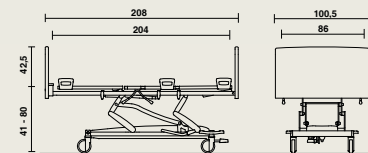

Materials and finishing

- Wood: Beech (massive, veneer), chipboard (melaminated).
- Metal: epoxy-coating RAL 9006 white aluminium.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials and construction**: please consult our material lists. We would refer to our extensive **colour gamut**.

Options

- 4-section mattress platform.
- Removable synthetic mattress platform sections.
- Foot extension up to + 20 cm.
- Panel types: Hypnos, Felina or Valentino.
- Bumpers (not possible on Midea).
- Wall spacer.
- Sideguards.
- (Reverse) Trendelenburg pedal.
- Cover cap.
- Wooden surrounds of the mattress platform frame.
- Decorative strip (front and back profile of the carriage).
- Accessory brackets.
- Integral S castors \varnothing 125 mm or \varnothing 150 mm on precision ball bearings, central brake.
- Hidden twin castors \varnothing 50 mm, central brake.
- Extendible drawer.

(*) Option

