

# Vico Hospital

**08700** Thanks to its hydraulic high/low adjustment, the Vico Hospital bed perfectly combines functionality, durability and an attractive price.

## Carriage

- Hydraulic high/low adjustment by means of supporting arms.
- The steel frame is finished with a removable vacuum-formed synthetic cover cap (\*).
- Trendelenburg adjustment up to  $\pm 12^\circ$  (\*).
- Reverse Trendelenburg adjustment up to  $\pm 13^\circ$  (\*).
- High-quality Levina design castors  $\varnothing$  125 mm, with separate brakes.
- Double design pedal to activate the central brake mechanism (\*).

## Mattress platform

- Three- or four- (\*) section mattress platform.
- The backrest can be adjusted up to  $\pm 70^\circ$ , with auto-regression by means of a gas spring.
- The knee rest can be adjusted up to  $\pm 34^\circ$  by means of a gas spring (in case of 4 sections).
- The legrest or footrest can be manually adjusted.
- The mattress platform is made of rounded steel mattress platform profiles.
- Standard equipped with cutouts for patients restraint straps.
- The backrest and footrest are equipped with 2 synthetic mattress holders.
- Standard with a telescopic bed extension up to + 20 cm.
- The 4 bed corners are equipped with an integrated bumper.
- The bed panels are standard equipped with a locking mechanism.
- The 4 bed corners are equipped with cases to install an IV rod or a lifting pole.
- Choice out of 3 different bed panels with a solid laminate filling panel.
- The side beams on both sides are equipped with 3 synthetic accessory hooks.

## Bed manipulations

- The high/low adjustment is activated by means of easily accessible pedals on both sides of the bed.
- The Trendelenburg position is activated by means of the foot pedal at the foot end (\*).
- The backrest and knee rest adjustments are supported by a gas spring. A lever on both bed sides guarantees a fluent activation.

### Paneltypes

- Vela
- Lauro (\*)
- Jaro (\*only with Trix)


### Sideguards

- Atmosphere (\*)
- Trix (\* only with Jaro-panel)

### Accessories

We would like to refer to our accessory list.


(\*) Option

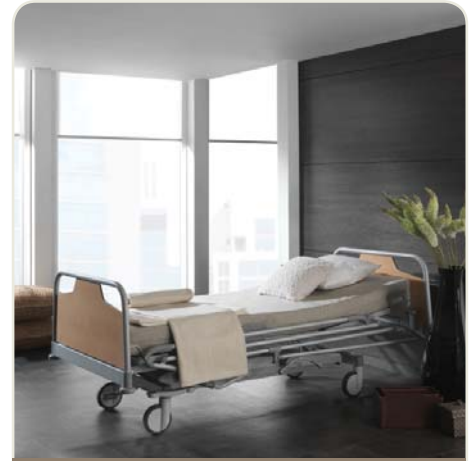

## Dimensions and technical data

- Width: 100 cm, Length: 212 cm
- Height variation of the mattress platform<sup>1</sup>:  $\pm 41$  cm -  $\pm 80$  cm
- Safe working weight: 200 kg
- Weight (excl. accessories): 115 kg
- Mattress platform dimensions: 204 x 86 cm
- Mattress dimensions: 195 x 85 cm, thickness: 12 à 16 cm

<sup>1</sup> With standard model

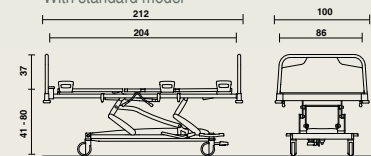

## Materials and finishing

- Metal: epoxy-coating RAL 9006 *white aluminium*, chrome (\*).
- Solid laminate.
- Self-lubricating synthetic hinges (= maintenance-free).
- Withstands the common cleaning products.

For more information about the **materials and construction**: please consult our material lists. We would refer to our extensive **colour gamut**.

## Options

- 4-section steel mattress platform.
- 3- or 4- section mattress platform with removable solid laminate sections of 6 mm.
- Foot extension up to + 20 cm.
- Panel types: Lauro or Jaro.
- Sideguards.
- (Reverse) Trendelenburg pedal.
- Removable vacuum-formed synthetic cover cap.
- Integral S castors  $\varnothing$  125 mm or  $\varnothing$  150 mm on precision ball bearings, centrally braked.
- Extendible drawer

